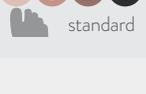
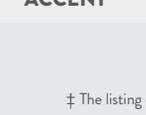
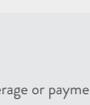


# College Park's Collection of Quality K3/K4 Feet

FOOT	MOUNTING	WEIGHT LIMIT	SIZES	IMPACT	HEIGHT	WEIGHT**	WARRANTY	L-CODE†	ENVIROHELL
 <b>SOLEUS</b>	Endo (4-hole)	220 lbs (21-24 cm) 250 lbs (25-26 cm) 275 lbs (27-30 cm)	21-30 cm		6.4-7.6 in	667 g	3	L5981	 narrow & wide
 <b>VELOCITY</b>	Endo (IP)	250 lbs (25-26 cm) 275 lbs (27-30 cm)	25-30 cm		4.7-5.2 in	631 g	3	L5981 + L5986	 narrow & wide
 <b>VENTURE</b>	Endo (IP) Exo (ALX)	250 lbs (21-24 cm) 275 lbs (25-30 cm)*	21-30 cm		2.5-2.7 in (IP) 5.0-5.2 in (ALX)	585 g	3	L5981 + L5986	 narrow & wide
 <b>TRUSTEP</b>	Endo (IP) Exo (ALX)	300 lbs (Endo) 352 lbs (Exo)	22-31 cm		3.5 in (IP) 6.0 in (ALX)	678 g	3 18 mo / 300+ lbs	L5979	 narrow & wide
 <b>HORIZON HD</b>	Endo (IP)	330 lbs (21-25 cm) 365 lbs (26-30 cm)	21-20 cm		2.3 in	559 g	3	L5981	 standard
 <b>ODYSSEY<sub>K3</sub></b>	Endo (IP)	275 lbs	21-30 cm		2.95-3.04 in	749 g	3	L5981 + L5968	 standard
 <b>TRIBUTE</b>	Endo (IP)	220 lbs	21-30 cm		2.4 in	529 g	2	L5979	 standard narrow & wide
 <b>ACCENT</b>	Endo (IP) D. Pylon (DP)	200 lbs (21-24 cm) 250 lbs (25-28 cm)	21-28 cm		3.3-3.9 in (IP) 7.8-8.5 (DP)	608 g (IP) 761 g (DP)	2	L5981 + L5986 + L5990	 narrow & wide

\* Low impact only

\*\*26 cm, including Enviroshell

Caucasian, Tan, Brown, Jet Black

† The listing of codes with these products should not be construed as a guarantee for coverage or payment. Ultimate responsibility for the coding of services/products rests with the individual practitioner.

College Park designs and manufactures a full line of anatomically correct, custom built prosthetic solutions. College Park is ISO 13485 certified with an extensive focus on quality testing, and proudly sources and produces all products entirely in the USA. With precise engineering and quality manufacturing, College Park is committed to innovating human locomotion solutions for users all over the world.



**college park**

TECHNOLOGY for the HUMAN RACE

800.728.7950 | 586.294.7950 | www.college-park.com



493 LIT BRO TS 210412

©2021, College Park Industries, Inc. Truststep, and Enviroshell are registered trademarks of College Park Industries, Inc.



**truststep**<sup>®</sup>

*the ultimate toolbox.*



Original and unrivaled, the Trustep sets the standard for unsurpassed comfort and natural-feeling mobility. It was one of the first prosthetic feet to provide anatomically correct movement in all planes, allowing for nearly the same vertical motion, rotation, and stability as an anatomical foot. With a full range of options available, the versatile design allows for complete customization.



## VERTICAL COMPLIANCE

The offset axis points in the Trustep allow the foot to flex and absorb vertical forces just like an anatomical foot.

## VERTICAL DISPLACEMENT



## TRI-AXIAL TECHNOLOGY

True multi-axial movement occurs on all three planes. The superior plantar flexion, dorsiflexion, transverse rotation, and inversion/eversion make the Trustep one of the most anatomically correct prosthetic feet — allowing for ease of movement, excellent ground compliance, and comfort.



## EVOLUTION

Since its release in 1991, the Trustep has undergone significant changes as new technology emerges. Some of these improvements include aircraft-grade polymer bumpers for unbeatable durability, noise abatement features, and a redesigned lighter, stronger heel.

# ANATOMICAL LOCOMOTION

FROM HEEL STRIKE TO TOE-OFF.



On heel strike, the rear bumper compresses to bring the toe down to foot flat. This gives higher rates of ground contact on declines for firm footing.



During midstance, the front and rear bumpers compress simultaneously to provide cushion. Dual hinge points allow the foot to maintain ground contact throughout the gait cycle.



To ready the foot for a smooth toe-off, the Truststep's front bumper compresses to absorb impact and lift the heel. The rollover shape and contoured toes assist in forward propulsion.



## PATIENT TESTIMONIAL



“After I put my legs on in the morning, I don't take them off until I go to bed in the evening. I've been told that I have a great gait. As a longtime amputee, I've worn many legs, and it hasn't always been like that. I started wearing my Truststeps in 2001 and haven't looked back.”

– Reggie, Truststep user

See more at [college-park.com/testimonials](http://college-park.com/testimonials)

## TECHNICAL SPECIFICATIONS

### FUNCTION

Tri-Axial

### MOUNTING

Endo (IP) Exo (ALX)

### SIZES

22-31 cm

### WEIGHT LIMIT

Endo (IP) 300 lbs (136 kg)

Exo (ALX) 352 lbs (160 kg)

### BUILD HEIGHT

Endo (IP) 3.5 in (8.7 cm)

Exo (ALX) 6.0 in (15.2 cm)

### WEIGHT\*

678 g

### WARRANTY

3 years (18 mo. 300+ lbs)

### IMPACT LEVEL

K3/4

### L CODES\*\*\*

L5979

### ENVIROHELL®†

narrow & wide

† Caucasian, Tan, Brown or Jet Black. \*26 cm model, including Enviroshell. \*\*\*The listing of codes with these products should not be construed as a guarantee for coverage or payment. Ultimate responsibility for the coding of services/products rests with the individual practitioner.

