

College Park's Collection of Quality K3/K4 Feet


FOOT	MOUNTING	WEIGHT LIMIT	SIZES	IMPACT	HEIGHT	WEIGHT**	WARRANTY	L-CODE†	ENVIROSHELL
 SOLEUS	Endo (4-hole)	220 lbs (21-24 cm) 250 lbs (25-26 cm) 275 lbs (27-30 cm)	21-30 cm		6.4-7.6 in	667 g	3	L5981	 narrow & wide
 VELOCITY	Endo (IP)	250 lbs (25-26 cm) 275 lbs (27-30 cm)	25-30 cm		4.7-5.2 in	631 g	3	L5981 + L5986	 narrow & wide
 VENTURE	Endo (IP) Exo (ALX)	250 lbs (21-24 cm) 275 lbs (25-30 cm)*	21-30 cm		2.5-2.7 in (IP) 5.0-5.2 in (ALX)	585 g	3	L5981 + L5986	 narrow & wide
 TRUSTEP	Endo (IP) Exo (ALX)	300 lbs (Endo) 352 lbs (Exo)	22-31 cm		3.5 in (IP) 6.0 in (ALX)	678 g	3 18 mo / 300+ lbs		 narrow & wide
 HORIZON HD	Endo (IP)	330 lbs (21-25 cm) 365 lbs (26-30 cm)	21-20 cm		2.3 in	559 g	3	L5981	 standard
 ODYSSEY_{K3}	Endo (IP)	275 lbs	21-30 cm		2.95-3.04 in	749 g	3	L5981 + L5968	 standard
 TRIBUTE	Endo (IP)	220 lbs	21-30 cm		2.4 in	529 g	2	L5979	 standard narrow & wide
 ACCENT	Endo (IP) D. Pylon (DP)	200 lbs (21-24 cm) 250 lbs (25-28 cm)	21-28 cm		3.3-3.9 in (IP) 7.8-8.5 (DP)	608 g (IP) 761 g (DP)	2	L5981 + L5986 + L5990 (IP) L5980 + L5986 + L5990 (DP)	 narrow & wide

* Low impact only

**26 cm, including Enviroshell



Caucasian, Tan, Brown, Jet Black

College Park designs and manufactures a full line of anatomically correct, custom built prosthetic solutions. College Park is ISO 13485 certified with an extensive focus on quality testing, and proudly sources and produces all products entirely in the USA. With precise engineering and quality manufacturing, College Park is committed to innovating human locomotion solutions for users all over the world.

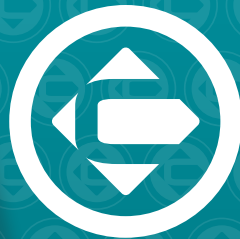




The first to incorporate Integrated Spring Technology (IST[®]), the Soleus provides natural movement with a progressive, 3-staged response and superior range of motion. The foot integrates Intelliweave composite springs through a tough, flexible plantar belt that responds to varying impacts. The unique, scaled design also provides 10 mm of vertical compliance for unbeatable comfort. The Soleus accommodates a moderate activity lifestyle up to the highest impact Paralympic athlete.

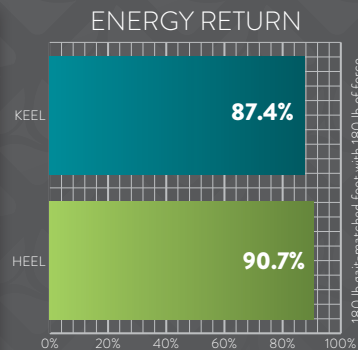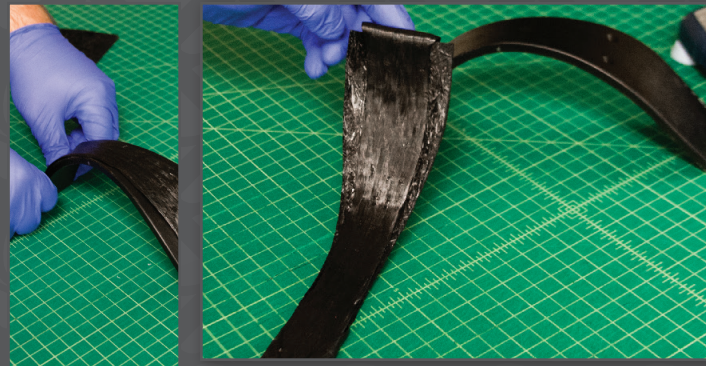

INTELLIWEAVE[®] COMPOSITE

College Park's proprietary Intelliweave composite technology brings the maximum amount of strength and flexibility. The fibers are woven by hand in a 3D pattern for ultimate durability and precision gait matching.

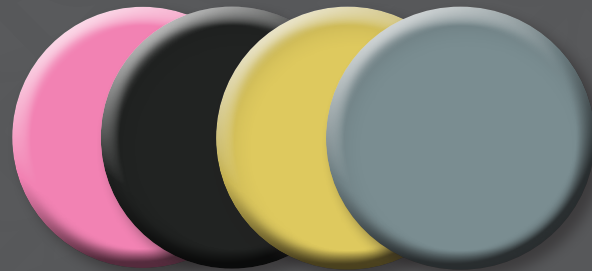

PERSONALIZED COLOR AND ENGRAVING

The Soleus housing is available in Perfect Pink, Midnight Black, Gold, and Silver and can also be engraved with one of nine different engravings to truly individualize the prosthesis.




WATERPROOF

Saltwater tested for 2 million cycles with an antimicrobial and UV inhibited foot shell, the progressive Soleus foot performs flawlessly at the pool, lake, ocean or outback.




INTEGRATED SPRING TECHNOLOGY[®]

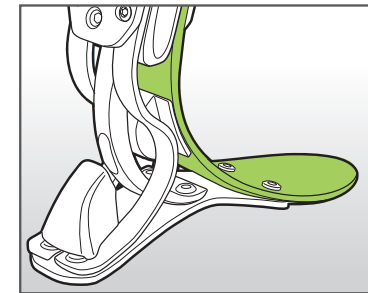


WALKING

The toe spring provides smooth movement and comfort for normal walking.

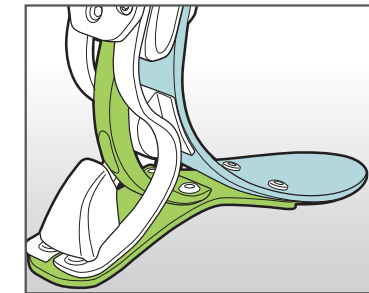

MODERATE ACTIVITY

During moderate activity, the midspring engages for a firmer response.

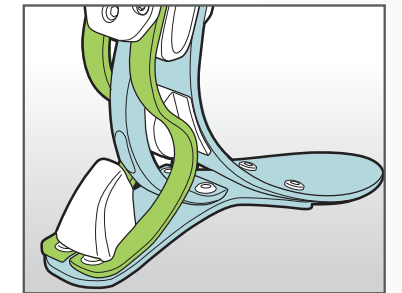

HIGH ACTIVITY

At the highest impact activity, the heel springs kick into action, for increased resistance and energy return.



"A lot of my military buddies use the Soleus. The Soleus is lightweight and gives me the great energy return I need. Plus it can be customized to fit with my style. It's my all-around everyday foot now, and I have the advantage of taking it with me into the water for watersports."
- Gideon, Soleus user. See more at www.college-park.com

TECHNICAL SPECIFICATIONS

FUNCTION

Integrated Spring Technology

MOUNTING

Endo (4-hole)

SIZES

21-30 cm

WEIGHT LIMIT

21-24 cm	220 lbs (100 kg)
25-26 cm	250 lbs (113 kg)
27-30 cm	275 lbs (125 kg)

CLEARANCE

21-24 cm	6.4 in (16.3 cm)
25-26 cm	7.1 in (18.0 cm)
27-30 cm	7.6 in (19.3 cm)

WEIGHT*

667 g

WARRANTY

3 years

L CODES**

L5981

IMPACT LEVEL



ENVIROSHELL[®]†



*26cm model, including Enviroshell. **The listing of codes with these products should not be construed as a guarantee for coverage or payment. Ultimate responsibility for the coding of services/products rests with the individual practitioner.
† Caucasian, Tan, Brown or Jet Black.

