

K-LEVEL
K3

CATEGORY
TRI-AXIAL



venture®

INTO THE
WILD




THE ULTIMATE FUSION




With a full length toe lever, the Venture provides a higher frequency dynamic response for more active users. Highly engineered bumpers and flexible Intelliweave® composites work together for ultimate progressive response. The tri-axial, custom built foot includes our exclusive Stride Control™ feature for effortless fine-tuning without disassembly.

2-STAGE RESPONSE

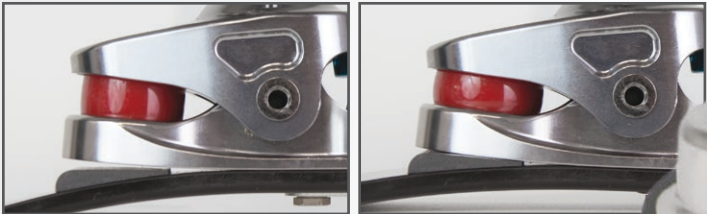

The front bumper provides comfort during walking, while the pocket limits the compression to engage the high-frequency composite keel during high activity.

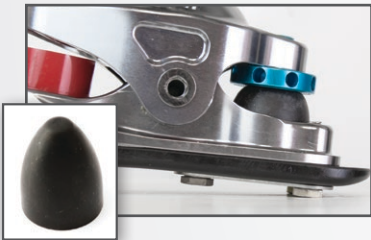

The rear bumper compresses at heel strike for cushioned comfort. The unique bullet shape provides controlled response for ground compliance.



TRI-AXIAL TECHNOLOGY

True multi-axial movement occurs on all three planes - the sagittal, coronal and transverse. The Venture ankle brings the most anatomical movement in a prosthetic foot, providing superior ground compliance.



CUSTOMIZABLE DESIGN

The Venture's bumpers and bushings can be swapped out as needed to increase and decrease resistance for a perfect feel. The bumpers were redesigned using aircraft-grade polymers for durability and longevity.



With only one pin to disassemble, the Venture requires minimal tools.



STRIDE CONTROL™

Stride Control allows for precise, incremental gait adjustments without any disassembly. Easily fine-tune the feel of the foot, increasing comfort and function for the user. A simple turn of a 4mm hex key changes the preload in the system for more or less resistance.



Adjusting the stride control up reduces heel resistance.

Adjusting Stride Control down will allow for a firmer heel response.

TECHNICAL SPECS

FUNCTION	MOUNTING	SIZES	IMPACT LVL	WEIGHT LIMIT	
Tri-Axial	Endo (IP) Exo (ALX)	21-30 cm		21-24 cm	250 lbs (113 kg)
				25-30cm	275 lbs** (125 kg)
CLEARANCE	WEIGHT*	WTY	L-CODE***	ENVIROSHELL†	
2.5-2.7 in (IP) (6.4-6.9 cm)	585g	3 years	L5981 L5986	 narrow & wide	
5.0-5.2 in (ALX) (12.7-13.2 cm)					

*26cm model, including Enviroshell. **Low impact only. ***The listing of codes with these products should not be construed as a guarantee for coverage or payment. Ultimate responsibility for the coding of services/products rests with the individual practitioner. † Caucasian, Tan, Brown or Jet Black.

ADDITIONAL ACCESSORIES



COLLEGE PARK SHELLTREAD
Provides additional grip and durability to Enviroshell.

VENTURE BUMPER KIT

Provides a toolbox of options for maintenance and to fine-tune the foot.



VENTURE

Read testimonials from users at www.IAmCollegePark.com





College Park designs and manufactures a full line of anatomically correct, custom built lower limb prosthetic solutions. College Park is ISO 13485 certified with an extensive focus on quality testing, and proudly sources and produces all products entirely in the USA. With precise engineering and quality manufacturing, College Park is committed to innovating human locomotion solutions for users all over the world.



college park
technology for the human race

800.728.7950 | 586.294.7950 | www.college-park.com



facebook.com/CollegeParkInd



youtube.com/CollegeParkInd



linkedin.com/company/College-Park-Industries